




058-732 BIODEX TABLE VASCULAR U/S VASC PRO 115V

058-732 Biodex Table Vascular
U/S Vasc PRO 115V
Item No: 058-732

The Vasc Pro Vascular Ultrasound Table is designed with ergonomics in mind.

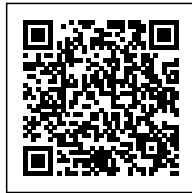

SKU: N/A

Price:

Stock: Available on backorder

Categories: [IMAGING TABLES](#), [Ultrasound Table](#)

Tags: [imaging table](#), [ultrasound table](#)

1

PRODUCT DESCRIPTION

058-732 Biodex Table Vascular
U/S Vasc PRO 115V
Item No: 058-732

The Vasc Pro Vascular Ultrasound Table is designed with ergonomics in mind. It addresses sonographer comfort and safety and proper patient positioning, virtually eliminating repetitive-strain injuries from vascular scanning.

Sonographer comfort and safety are accomplished through features such as the new adjustable foot support that accommodates patients of various heights, and optional side rails enabling close access to the patient without sonographer contortion.

Fowler positioning to 80 degrees, which accommodates natural body extension and sitting



position without slide or shift, an optional Articulating Scanning Arm Board adjustable from 0 to 130-degrees locking every 10 degrees, the thick antimicrobial mattress and table adjustability all create a comfortable imaging experience for the patient.

Suitable Procedures

Lower Extremity

Venous Mapping and Doppler (DVT)

Venous Reflux

Venous Insufficiency

Venous Ablation

Arterial Scan/Doppler (PAD)

Upper Extremity

Arterial*

Venous Access Placement and Assessment

Carotid

Transcranial

Aorta

Breast

Renal Transplant

Testicular

SMA

Carotid Duplex

Abdomen

Thyroid

Renal

Renal Artery

Musculoskeletal

*Requires optional vascular arm board.

All motorized tables are ETL product listed to the UL 60601-1 Medical Electrical Equipment Safety Standard.