



PATRIOT HOMECARE BEDS SEMI-ELECTRIC

Patriot Homecare Beds Semi-Electric

Total bed weight assembled: 166 lb (bed and motor only).

Available with Fiberboard Decorative Walnut Bed-ends or High-impact Plastic Bed-ends.


Maximum patient weight 350 lbs.

450 lb evenly distributed total weight capacity including patient, mattress, rails, and accessories.

Bed Height: Min. 15", Max. 24" (at deck level).

Lift Range: Head - 60°, Knee - 30°.

Overall Bed Size: 87" L x 36" W.

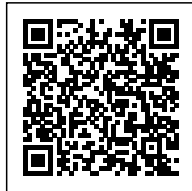

SKU: N/A

Price:

Stock: 10 in stock (can be backordered)

Categories: [BED](#), [Semi Electric Beds](#)

Tags: [Homecare bed](#), [patient bed](#), [semi electric bed](#)

ADDITIONAL INFORMATION

Item No. 2	1738, 1739, 1740, 1741, 1742, 1743, 1744, 1745, 1746, 1747, 1748, 1749, 1750, 1751, 1752, 1753, 1754, 1755, 1756, 1757, 1758, 1759, 1760
-------------------	---

VARIATIONS







Image	SKU	Price	Description	Item No. 2
				US0208
				US0208-220V
				US0208-220V
				US0208-RFPKG



Image	SKU	Price	Description	Item No. 2
				US0208-RFPKGHR
				US0208-RPKG
				US0208-RPKGHR
				US0208PL



Image	SKU	Price	Description	Item No. 2
				US0208PL-220V
				US0208PL-RFPKG
				US0208PL-RFPKGHR
				US0208PL-RPKG



Image	SKU	Price	Description	Item No. 2
				US0208PL-RPKGHR
				US0218-IFPKGFR
				US0218-IFPKGHR
				US0218-RFPKG



Image	SKU	Price	Description	Item No. 2
				US0218-RFPKGHR
				US0218-RPKG
				US0218-RPKGHR
				US0218PL-RPKG



Image SKU Price Description Item No. 2



US0218PL-RFPKGHR



US0218PL-RPKG



US0218PL-RPKGHR

PRODUCT DESCRIPTION

Patriot Homecare Beds Semi-Electric

Introducing the new Patriot Semi-Electric Homecare Bed with a grid sleep surface!
Designed and engineered by the same team who develop the GF Long Term Care beds.
Sleep surface and motors are MADE IN THE USA!

The motor unit and bed design provide the ultimate in quiet, smooth operation.
Complete welded steel frame construction.

The head and foot sections can be raised or lowered simultaneously or independently with one hand control.

During a power failure the motor unit can lower both the head and foot section with just a



single standard 9 volt battery.

Assembles easily with just five components.

Designed to meet the requirements of applicable HCPCS codes.

Total bed weight assembled: 166 lb (bed and motor only).

Available with Fiberboard Decorative Walnut Bed-ends or

High-impact Plastic Bed-ends.

Maximum patient weight 350 lbs.

450 lb evenly distributed total weight capacity including patient, mattress, rails, and accessories.

Bed Height: Min. 15", Max. 24" (at deck level).

Lift Range: Head - 60°, Knee - 30°.

Overall Bed Size: 87" L x 36" W.

No external pinch points!

Slots in frame to ensure proper rail placement.

Locking pins secure bed ends to sleep deck. Attached lanyards prevent loss while moving or storing.

Only homecare bed to offer a 4-inch extension kit that can extend either the head, foot, or both head and foot decking. Bed length will be 84" with one kit, 88" with 2 kits.

The only homecare bed to offer optional corner mattress retainers .

Grid Deck offers the following advantages:

Much stronger frame.

Better mattress support.

No spring sag or bottoming out with heavier patients.

Easier to clean.

No connecting links or springs to lose.

Warranty:

Two year warranty for defects in workmanship and materials of mechanical components, frame, and electronics.

During the warranty period, defective items will be repaired or replaced at manufacturer's option. The warranty does not include any labor charges incurred in replacement part(s) installation or any associated freight or shipping charges to the manufacturer.

Bed shipped in three cartons:

1. Foot section and motor.
2. Head section and Hi/Lo rod.
3. Bed end set.

US0208 Patriot Semi-Electric Bed Only

US0208-220V Patriot Semi-Electric Bed Only - 220 V - International Only



US0208-IFPKGHR Patriot Semi-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and Half Chrome Rails

US0208-RFPKG Patriot Semi-Electric Bed with 1633-Reversible Foam Mattress and Full Chrome Rails

US0208-RFPKGHR Patriot Semi-Electric Bed with 1633-Reversible Foam Mattress and Half Chrome Rails

US0208-RPKG Patriot Semi-Electric Bed with No Mattress and Full Chrome Rails

US0208-RPKGHR Patriot Semi-Electric Bed with No Mattress and Half Chrome Rails

US0208PL Patriot Semi-Electric Bed Only with Plastic Ends

US0208PL-220V Patriot Semi-Electric Bed Only - 220 V - International Only, Plastic Ends

US0208PL-RFPKG Patriot Semi-Electric Bed with 1633-Reversible Foam Mattress and Full Chrome Rails, Plastic Ends

US0208PL-RFPKGHR Patriot Semi-Electric Bed with 1633-Reversible Foam Mattress and Half Chrome Rails, Plastic Ends

US0208PL-RPKG Patriot Semi-Electric Bed with No Mattress and Full Chrome Rails, Plastic Ends

US0208PL-RPKGHR Patriot Semi-Electric Bed with No Mattress and Half Chrome Rails, Plastic Ends

US0218-IFPKGFR Patriot Semi-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and FDA Full Rails

US0218-IFPKGHR Patriot Semi-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and FDA Half Rails

US0218-RFPKG Patriot Semi-Electric Bed with 1633- Rolled Foam Mattress and FDA Full Rails

US0218-RFPKGHR Patriot Semi-Electric Bed with 1633- Rolled Foam Mattress and FDA Half Rails

US0218-RPKG Patriot Semi-Electric Bed with No Mattress and FDA Full Rails

US0218-RPKGHR Patriot Semi-Electric Bed with No Mattress and FDA Half Rails

US0218PL-RFPKG Semi-Electric Bed with 1633- Rolled Foam Mattress and FDA Full Rails, Plastic Ends

US0218PL-RFPKGHR Semi-Electric Bed with 1633- Rolled Foam Mattress and FDA Half Rails, Plastic Ends

US0218PL-RPKG Semi-Electric Bed with No Mattress and FDA Full Rails, Plastic Ends

US0218PL-RPKGHR Semi-Electric Bed with No Mattress and FDA Half Rails, Plastic Ends

ADDITIONAL INFORMATION



Item No. 2

1738, 1739, 1740, 1741, 1742, 1743, 1744, 1745, 1746, 1747,
1748, 1749, 1750, 1751, 1752, 1753, 1754, 1755, 1756, 1757,
1758, 1759, 1760