



## PATRIOT LX SEMI-ELECTRIC HOME-CARE BED

Patriot LX Semi-Electric Home-care Bed

Technical Specifications


Total bed weight assembled: 118 lb

Maximum patient weight: 350 lbs EVENLY DISTRIBUTED, 450 lb distributed total weight capacity including patient, mattress, rails and accessories

Overall Bed Size: 84" L x 35" W

Lift Range: Head - 60, Knee - 30

Bed Height: Min. 10", Max. 20" (at deck level)

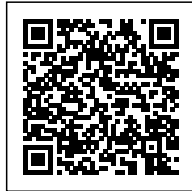

**SKU:** N/A

**Price:**

**Stock:** N/A

**Categories:** [BED](#), [Semi Electric Beds](#)

**Tags:** [Homecare bed](#), [patient bed](#), [semi electric bed](#)

### ADDITIONAL INFORMATION

**Item No. 2** 1709, 1710, 1711, 1768, 1712, 1767,  
1713, 1714, 1715, 1716, 1769, 1717

### VARIATIONS



Image	SKU	Price	Description	Item No. 2
				US5000
				US5000-RFPKG
				US5000-RFPKGHR
				US5000-RPKG



<b>Image</b>	<b>SKU</b>	<b>Price</b>	<b>Description</b>	<b>Item No. 2</b>
				US5000-RPKGHR
				US5000-RPKGQR
				US5000-SFPKG
				US5000-SFPKGHR



<b>Image</b>	<b>SKU</b>	<b>Price</b>	<b>Description</b>	<b>Item No. 2</b>
				US5000-RFPKGQR
				US5000-RFPKGC
				US5000-SFPKGC
				690-8084-000

**PRODUCT DESCRIPTION**

Patriot LX Semi-Electric Home-care Bed

Robust steel grid deck sleep surface provides superior mattress support and easy to clean



surface compared to traditional spring decks

Head and foot sections can be raised or lowered simultaneously or independently with one hand control

Recessed main frame offers maximum pinch point protection between grid deck and frame - reducing the risk of injury

Lightweight design weighs only 118 lbs - lightest bed on the market

Bed setup is so simple that it can save on setup time and requires minimal maintenance keeping your technicians time efficient

Patent pending integrated HI-LO mechanism makes bed height adjustment easy and reduces potential loss of parts such as HI/LO rod when used as rentals

Self contained crank for semi-electric that can be used at either end makes adjustment from low 10" to 20" easy and even allows for trendelenburg and slant positioning

Patent pending rail slots in frame ensure quick and proper rail placement

Single motor unit is designed for easy installation, requires less maintenance and provides the ultimate in quiet, smooth operation

Attractive, "home-like" thermo fused bed ends offer a less institutional look and eliminate a bed component that often requires frequent repair and maintenance

Features SilverSolutions- an anti-microbial technology built into the paint that protects against bacteria and odor

During a power failure, the motor can lower both the head and foot section with just a single 9 volt battery

4" extension kit can be installed on head, foot or both to extend bed length to 84" or 88"

Grid Deck offers the following advantages:

Much stronger frame

Better mattress support

No spring sag or bottoming out with heavier patients

Easier to clean

No connecting links or springs to lose

Technical Specifications

Total bed weight assembled: 118 lb

Maximum patient weight: 350 lbs EVENLY DISTRIBUTED, 450 lb distributed total weight capacity including patient, mattress, rails and accessories

Overall Bed Size: 84" L x 35" W

Lift Range: Head - 60, Knee - 30

Bed Height: Min. 10", Max. 20" (at deck level)

US5000 Patriot LX Semi-Electric Bed Only



US5000-RFPKG Patriot LX Semi-Electric Bed with 1633-Rolled Foam Mattress and FDA Full Rails

US5000-RFPKGHR Patriot LX Semi-Electric Bed with 1633-Rolled Foam Mattress and FDA Half Rails

US5000-RFPKGQR Patriot LX Semi-Electric Bed with 1633-Rolled Foam Mattress and FDA Quarter Rails

US5000-RFPKGC Patriot LX Semi-Electric Bed with 1633-Rolled Foam Mattress and Chrome Full Rails

US5000-RPKG Patriot LX Semi-Electric Bed with No Mattress and FDA Full Rails

US5000-RPKGHR Patriot LX Semi-Electric Bed with No Mattress and FDA Half Rails

US5000-RPKGQR Patriot LX Semi-Electric Bed with No Mattress and FDA Quarter Rails

US5000-SFPKG Patriot LX Semi-Electric Bed with 1633-Standard Foam Mattress & FDA Full Rails

US5000-SFPKGHR Patriot LX Semi-Electric Bed with 1633-Standard Foam Mattress & FDA Half Rails

US5000-SFPKGC Patriot LX Semi-Electric Bed with 1633-Standard Foam Mattress & Chrome Full Rails

690-8084-000 Extension Kit for Patriot LX

## **ADDITIONAL INFORMATION**

**Item No. 2** 1709, 1710, 1711, 1768, 1712, 1767, 1713, 1714, 1715, 1716, 1769, 1717