



UNIVERSAL SLINGS

Universal Slings



MS882 SURELIFT™ Universal Mesh Sling - Medium

MS883 SURELIFT™ Universal Mesh Sling - Large

UP892 SURELIFT™ Universal Sling - Medium

UP893 SURELIFT™ Universal Sling - Large

UP894 SURELIFT™ Universal Sling - XX-Large, Bariatric

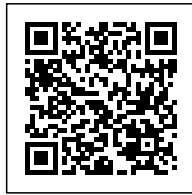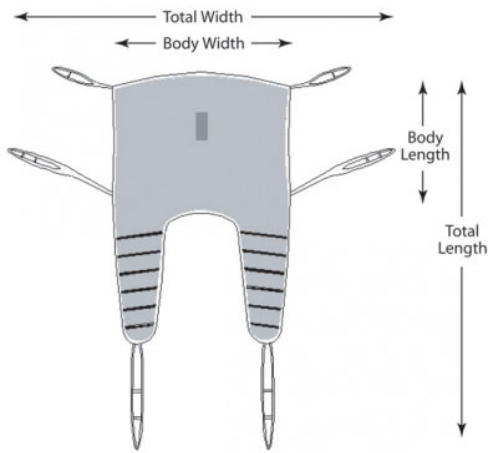

SKU: N/A

Price:

Stock: 10 in stock (can be backordered)

Categories: [Lift Slings](#), [PATIENT LIFTS](#)

Tags: [lift slings](#), [slings](#)

ADDITIONAL INFORMATION

Weight N/A

Dimensions N/A

Item No. 2 392, 393, 394, 395, 396

VARIATIONS




Image

SKU


Price

Description


Item No. 2

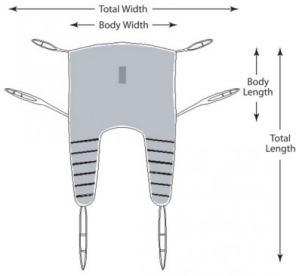


MS882

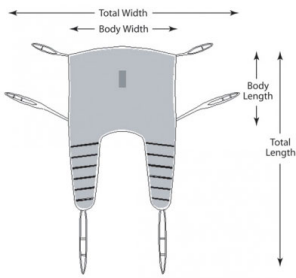

MS883

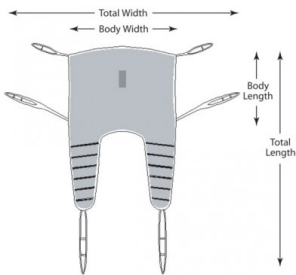

UP892

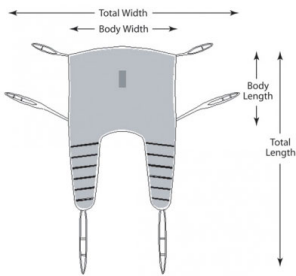

UP893




Image

SKU Price

Description

Item No. 2

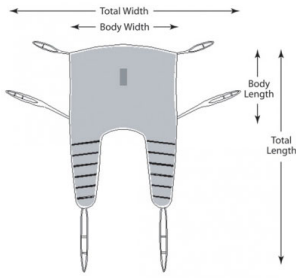

UP894

PRODUCT DESCRIPTION

Universal Slings

MS882 SURELIFT™ Universal Mesh Sling - Medium

MS883 SURELIFT™ Universal Mesh Sling - Large

UP892 SURELIFT™ Universal Sling - Medium

UP893 SURELIFT™ Universal Sling - Large

UP894 SURELIFT™ Universal Sling - XX-Large, Bariatric

Universal slings (also called divided leg slings) are the easiest sling to use while transferring a resident from a seated position. These can also be used for lifting from a supine position, making them suitable for most lifting situations. They are an easy-to-fit general purpose sling designed for the majority of users. They can be used for transfers from a bed to a chair, wheelchair, or toilet, chair to chair, or from the floor to a bed or chair.

Universal slings should be used for residents that are totally or partially dependent, non-weight bearing, or bariatric. Those who have limited head control should use a universal sling with head support. The large opening accommodates toileting. Mesh slings should be used for bathing. LUMEX Universal slings have four or six connection points and can be used with 2-, 4-, or 6-point spreader bars.

- Provides excellent patient comfort and security
- Reinforced back for support
- Non-mesh slings have padding on legs and back
- Height-Adjustable straps
- Easy to use
- Machine washable



- Six point connectivity is designed to enhance positioning and safety
- Designed to meet the requirements of HCPCS Code: E0621
- Six (6) month limited warranty

ADDITIONAL INFORMATION

Weight	N/A
Dimensions	N/A
Item No. 2	392, 393, 394, 395, 396